

Prosthetic components **Torque**

ANKYLOS® Cover Screw



manual

ANKYLOS® Cover Screw for SynCone® C/
Abutment straight 5°



manual

ANKYLOS® Gingiva Former



manual

ANKYLOS® Balance C/
Temporary Abutment



manual

Prosthetic components **Torque**

ANKYLOS® Balance Anterior C/
and /X Abutment



15 Ncm

ANKYLOS® Balance Posterior C/ Abutment



15 Ncm

ANKYLOS® CERCON® Balance C/ Abutment




15 Ncm

ANKYLOS® Regular C/ and /X Abutment



15 Ncm

ANKYLOS® TitaniumBase C/ and /X
ANKYLOS® ScanBase C/ and /X

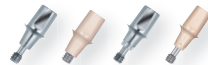

15 Ncm

ANKYLOS® SynCone® C/ Abutment



15 Ncm

ATLANTIS™ Abutment for ANKYLOS®



15 Ncm

Prosthetic components **Torque**

ANKYLOS® Abutment with
Snap Attachment C/



25 Ncm

ANKYLOS® LOCATOR® Abutment C/



30 Ncm

Prosthetic components **Torque**

ANKYLOS® Balance Base Abutment C/
(SmartFix™)



straight 25 Ncm
angled 15 Ncm (abutment body),
angled 25 Ncm (abutment head)

ANKYLOS® Standard C/ Abutment



straight 25 Ncm
angled 15 Ncm



ANKYLOS® Hexagon Screwdriver
1.0 mm
Order no. 3103 3400



ANKYLOS® Insert for Prosthetic Ratchet
Order no. 3103 3625
1.0 mm Hex, 15 Ncm
Order no. 3103 3628 (not shown)
1.0 mm Hex, long, 15 Ncm



ANKYLOS® Insert for Prosthetic Ratchet for LOCATOR®
30 Ncm
Order no. 3103 3635



1


ANKYLOS® Insert for Prosthetic Ratchet for Balance Base Abutment

1 | 1.0 mm hex, 15 Ncm
Order no. 3103 3625 short
Order no. 3103 3628 long (not shown)



2

2 | 1.8 mm hex, 25 Ncm
Order no. 3103 3627

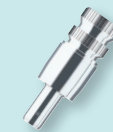

1

ANKYLOS® Insert for Prosthetic Ratchet for Standard Abutment

Straight
1 | for a-Abutment, 25 Ncm
Order no. 3103 3621 short
Order no. 3103 3633 long (not shown)



2

Straight
2 | for b-Abutment, 25 Ncm
Order no. 3103 3622 short
Order no. 3103 3634 long (not shown)



3

Angled
3 | 1.0mm Hex, 15 Ncm
Order no. 3103 3625 short
Order no. 3103 3628 long (not shown)